

CAMouflage[®] NOW

“The Speed You Desire”



The combination of **speed**, **esthetics** and **durability** you desire

- Indicated for inlays, onlays, and single-unit anterior and posterior crowns
- Mill in under 10 minutes*
- No firing, staining or glazing needed prior to delivery
- Easy finish, with high gloss retention over time
- Exclusively available for Glidewell milling systems



RETAINS
SHADE



STAIN
RESISTANT

Before



After

This patient's onlay on tooth #19 and inlay on tooth #20 were milled in minutes from CAMouflage NOW.

“CAMouflage[®] NOW is an excellent restorative alternative for indirect dentistry. I can provide quality, conservative restorations to my patients in a matter of minutes. This, in turn, has increased the productivity of my practice.”

– Chad Duplantis, DDS | Fort Worth, TX

For more information
glidewelldental.com ■ 800-854-7256



**GLIDEWELL
DENTAL**


BENEFITS

- Rapid patient delivery with no oven wait
- Integrates easily into your workflows
- Shade stability long after delivery
- Natural opalescence and fluorescence for maximum esthetics
- Strong and long-lasting restorations, with enhanced resistance to chipping at the margin
- Affordable pricing, billable to insurance as ceramic



Sold in 5-block packs

READY TO DELIVER IN UNDER 15 MINUTES

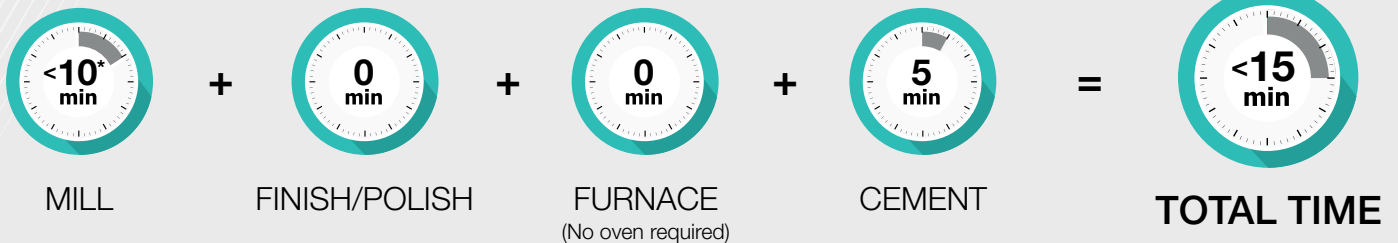

SPECIFICATIONS

MATERIAL

- Approximately 80 percent ceramic-filled nanohybrid composite
- Flexural Strength: 193 MPa

INDICATIONS


- Inlays, onlays, anterior and posterior crowns

SHADES

- Available in seven shades that correspond to the VITA® classical shade system, five of which are available in two translucencies, and one bleach shade: A1HT, A1LT, A2HT, A2LT, A3HT, A3LT, A3.5HT, A3.5LT, B1HT, B1LT, C2LT, D2LT, BLLT

EASY PREPARATION

ONLY

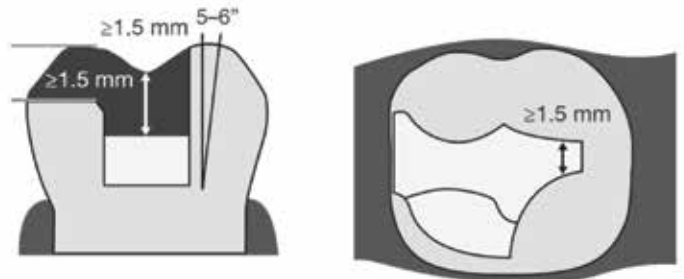

- Prepare cavity walls with tapered 5–6° to the long axis.
- Wall thickness should be ≥ 1.5 mm.
- Isthmus thickness should be ≥ 1.5 mm.

*Milling times vary based on prescription and mill used. Average milling time is shown.

VITA is a registered trademark of VITA Zahnfabrik.

MKT 10487_2.0 PK-280353-092018

For more information
glidewelldental.com ■ 800-854-7256

